

# **Cooling incubators**







**KB**| **KT** series

# Cooling incubators

Series KT	Page 6
Series KB	Page 8
 Options and Accessories	Page 11

BINDER cooling incubators are gentle on samples and provide high-quality performance. Every KB Series cooling incubator thus makes safe, reproducible incubation possible, even at high ambient temperatures. Together with the efficient weekly program function, a cooling incubator from BINDER is truly versatile in laboratory work with microorganisms. The KT Series combines exceptional performance with low energy consumption. Additionally, every BINDER cooling incubator offers true added value with its exemplary environmental sustainability.

# PRODUCT SELECTION COOLING INCUBATORS

Features	KT > Page 6	KB ⇒Page 8
Model		
Model 23	-	•
Model 53	•	•
Model 115	•	•
Model 170	•	-
Model 240	-	•
Model 400	-	•
Model 720	_	•
General data		
Temperature range [°C]	4100	-5100
Ethernet interface	•	•
USB interface	•	•
Cooling technology	Thermoelectric cooling	Compressor
Controller and timer functions		
Controller display	LCD	
Programming option	•	•
Weekly program function	•	•
Adjustable fan speed	•	•
Real-time clock	•	•
Class 1 temperature safety device	•	•
Class 3.1 independent temperature safety device	•	•
Temperature alarm, acoustic	•	•
Temperature alarm, visual	•	•
•Standard equipment		

•Standard equipment

## ADVANTAGES AT A GLANCE



HOMOGENEOUS TEMPERATURE DISTRIBUTION

#### Series KT, KB


BINDER APT.line<sup>™</sup> preheating chamber technology ensures long-term operational stability even under full load.



#### EASY DECONTAMINATION

Series KT, KB

The disinfection routine at 100 °C and high quality stainless-steel inner chamber makes it easy to clean.

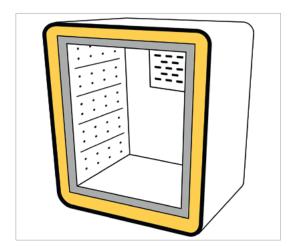

#### EXTREMELY ROBUST AND LONG-LASTING QUALITY

#### Series KT, KB

HIT insulation provides a layer of foam between the boiler and the housing. This guarantees maximum stability and particle-free thermal insulation with no thermal bridges.



FLEXIBLE, MIXED STACKING OPTIONS

Series KT

Both units have identical footprints in sizes KT 115 and KT 170.

# Series KT | with thermoelectric cooling

The KT series combines outstanding performance with impressive energy efficiency and environmental friendliness. The cooled incubators of the KT series are state-of-the-art in application and protect the samples.





KT 115 model

#### MAIN FEATURES

- Temperature range: 4 °C to 100 °C
- APT.line<sup>™</sup> preheating chamber technology
- Adjustable fan speed
- Electric cooling with Peltier module
- Controller with time-segment and real-time programming

#### AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
KT 53	53	› go2binder.com/en-KT53
KT 115	102	› go2binder.com/en-KT115
KT 170	163	> go2binder.com/en-KT170

- Display via LCD monitor
- Input via pushbutton/rotary knob
- Inner door made of tempered safety glass
- 2 stainless steel racks
- Units up to 115 liters are stackable

#### TIPS & TRICKS

i

Benefit from our many years of experience, as demonstrated in our technical papers:

Tips & Tricks: Case Studies: > go2binder.com/en-tips-tricks
> go2binder.com/en-case-studies



Find all details online go2binder.com/en-KT

#### **TECHNICAL DATA**

Description	KT 53	KT 115	KT 170
Measures - Outer dimensions			
Width net [mm]	660	860	860
Height net [mm]	635	715	1025
Depth net [mm]	630	655	655
Measures - Internal Dimensions			
Width [mm]	400	600	600
Height [mm]	400	455	765
Depth [mm]	334	355	355
Measures			
Interior volume [L]	53	102	163
Net weight of the unit (empty) [kg]	63	83	102
permitted load [kg]	40	100	120
Load per rack [kg]	15	30	30
Performance Data Temperature			
Temperature range [°C]	4100	4100	4100
Temperature variation at 37 °C [± K]	0.3	0.3	0.4
Temperature fluctuation [± K]	0.1	0.1	0.1
Recovery time after 30 seconds door open at 37 °C [min]	2	3	3
Electrical data			
Nominal power (version 120 V) [kW]	0.4	0.7	0.8
Nominal power (version 230 V) [kW]	0.4	0.7	0.8
Fixtures			
Number of shelves (std./max.)	2/5	2/6	2/10


#### ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	KT 53	KT 115	KT 170
200240 V 1 50/60 Hz	Standard	9020-0311	9020-0313	9020-0289
100120 V 1 50/60 Hz	Standard	9020-0312	9020-0314	9020-0310

▶ You will find a selction of options, accessories and services suitable for each product on page 11

#### DIMENSIONS incl. fittings and connections [mm]

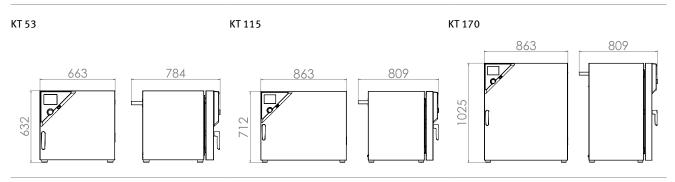

# Series KB | with compressor technology

The most versatile among the cooled incubators for microorganisms: The BINDER KB series cooled incubator controls temperature ranges of -5 °C to 100 °C. With its comprehensive program functions, this cooled incubator offers a wide range of capabilities and delivers reproducible test results.



KB 53 model

#### MAIN FEATURES

- Temperature range (Model 23): 0 °C to 100 °C
- Temperature range (Models 53, 115, 240, 400, 720):
   -5 °C to 100 °C
- APT.line<sup>™</sup> preheating chamber technology
- Cooling with compressor cooling unit
- Door heating system from 240 liters

#### AVAILABLE SIZES

Model	Interior volume [L]	Online data sheets
KB 23	20	› go2binder.com/en-KB23
KB 53	53	› go2binder.com/en-KB53
KB 115	115	› go2binder.com/en-KB115
KB 240	247	› go2binder.com/en-KB240
KB 400	400	> go2binder.com/en-KB400
KB 720	698	› go2binder.com/en-KB720



- Adjustable fan speed
- Controller with time-segment and real-time programming
- Inner door made of tempered safety glass
- 2 stainless steel racks
- Units up to 115 liters are stackable

#### TIPS & TRICKS

Benefit from our many years of experience, as demonstrated in our technical papers:

Tips & Tricks:	>go2binder.com/en-tips-tricks
Case Studies:	>go2binder.com/en-case-studies

Find all details online go2binder.com/en-KB

#### **TECHNICAL DATA**

Description	KB 23	KB 53	KB 115	KB 240	KB 400	KB 720
Measures - Outer dimensions						
Width net [mm]	435	635	835	925	925	1250
Height net [mm]	620	835	1025	1465	1950	1925
Depth net [mm]	520	580	650	800	805	885
Measures - Internal Dimensions						
Width [mm]	222	400	600	650	650	970
Height [mm]	330	400	480	785	1270	1250
Depth [mm]	277	330	400	485	485	576
Measures						
Interior volume [L]	20	53	115	247	400	698
Net weight of the unit (empty) [kg]	44	78	106	171	221	304
permitted load [kg]	25	40	50	100	120	150
Load per rack [kg]	10	15	20	30	30	45
Performance Data Temperature						
Temperature range [°C]	0100	-5100	-5100	-5100	-5100	-5100
Temperature variation at 37 °C [± K]	0.3	0.2	0.2	0.2	0.2	0.2
Temperature fluctuation at 37 °C [± K]	0.1	0.1	0.1	0.1	0.1	0.1
Recovery time after 30 seconds door open at 37 °C [min]	2	2	2	2	4	3
Electrical data						
Nominal power (version 120 V) [kW]	0.3	0.6	0.7	1.2	1.4	-
Nominal power (version 230 V) [kW]	0.3	0.6	0.7	1.2	1.4	2.3
Nominal power (version 240 V) [kW]	-	-	-	-	-	2.3
Fixtures						
Number of shelves (std./max.)	2/3	2/4	2/5	2/9	2/15	2/15


#### ORDERING INFORMATION [Art. No.]

Further details: Go to www.binder-world.com > Search > Enter Article Number

Voltage	Option model	KB 23	KB 53	KB 115	KB 240	KB 400	KB 720
230 V 1 50/60 Hz	Standard	9020-0112	9020-0199	9020-0397	-	-	-
200240 V 1 50 Hz	Standard	-	-	-	9020-0202	9020-0203	9020-0204 9020-0306
120 V 1 60 Hz	Standard	9020-0113	9020-0302	9020-0303	-	_	-
100120 V 1 60 Hz	Standard	_	_	-	9020-0304	9020-0305	_

 $\blacktriangleright$  You will find a selction of options, accessories and services suitable for each product on page 11

#### DIMENSIONS incl. fittings and connections [mm]

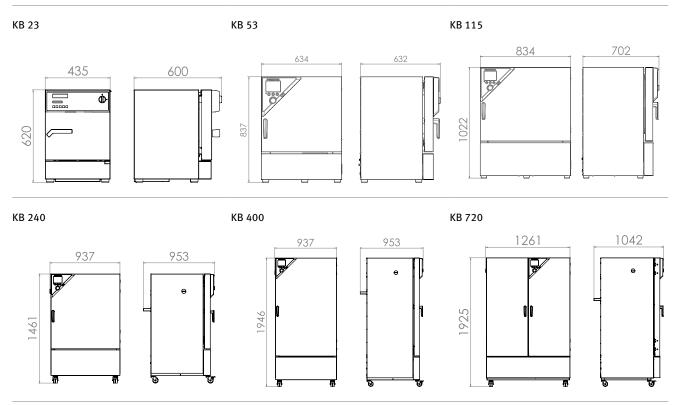

# **EXCERPT FROM OPTIONS AND ACCESSORIES**

Access port with silicone plug	For introducing external measuring instruments. Large selection of different diameters and posi- tions (left, right, or top) Series KT, KB
Interior socket	Socket with cover, with associated LEMO connector, nominal voltage, max. 500 W, max. 90 °C, and protection class IP67. For installing external units in the interior. Series KT, KB
Door lock	Lockable door handle prevents unauthorized access. Series KT, KB
Data Logger Kit T	Independent temperature data logging. The tailored product solution contains helpful accessories: From the magnetic mounting of the data logger to the BINDER unit and securing sensor cable for accurate readings. Series KT, KB
Services	You can access our installation, maintenance, calibration, validation, or warranty services by taking out a Maintenance Service Agreement or by simply getting in touch when you need us. More infor- mation is available in the accessories and services chapter. Series KT, KB
BINDER INDIVIDUAL	Using a BINDER series production product as a basis, we develop individual solutions as a single unit or in small batches. All products are tested, certified, and supplied with a full functional gu- arantee and warranty. More information is available in the accessories and services chapter. Series KT, KB









Access port with silicone plug

Shaker operation in the chamber

Door lock

Data Logger for temperature

### **BINDER** products

CO₂ incubato	rs
Incubators	
Cooling incut	
Growth cham	bers
Ultralow tem	perature freezers
Drying and h	eating chambers
Safety drying	g chambers
Vacuum dryir	ıg chambers
Constant clim	nate chambers
Dynamic clim	ate chambers
	Never miss a product novelty - Subscribe and keep up to speed with the BINDER monthly newsletter.
∎r <del>2</del> 59	Subscription
go2binder.co	m/en-newsletter-subscription



#### Best conditions for your success

#### BINDER GmbH

Tuttlingen, Germany

Tel +49 7462 2005 0 Fax +49 7462 2005 100 info@binder-world.com www.binder-world.com BINDER Asia Pacific (Hong Kong) Ltd. Kowloon, Hong Kong, P.R. China

Tel +852 39070500 Fax +852 39070507 asia@binder-world.com www.binder-world.com

BINDER Environmental Testing Equipment (Shanghai) Co., Ltd. Shanghai, P.R.China

Tel +86 21 685 808 25 Fax +86 21 685 808 29 china@binder-world.com www.binder-world.com Representative Office for CIS Moscow, Russia

Tel +7 495 988 15 16 Fax +7 495 988 15 17 russia@binder-world.com www.binder-world.ru

BINDER Inc. Bohemia, NY, USA

Tel +1 631 224 4340 Fax +1 631 224 4354 usa@binder-world.com www.binder-world.us